


abberior FLUX 660, SNAP, 30 nmol

Item number

FX660-0145-30NMOL

Description

abberior FLUX 660 is a dye developed for MINFLUX imaging in the red spectral region. It combines exceptional photostability, high brightness, and low background labeling with photoswitching kinetics being optimized for MINFLUX imaging. abberior FLUX 660 can be applied with excitation laser wavelengths between 630 - 670 nm together with a detection window between 650 - 750 nm. For MINFLUX imaging abberior FLUX 660 is used with MINFLUX imaging buffers such as GLOX buffer. abberior FLUX 660 can substitute dyes like Alexa Fluor® 660, or CF660®.

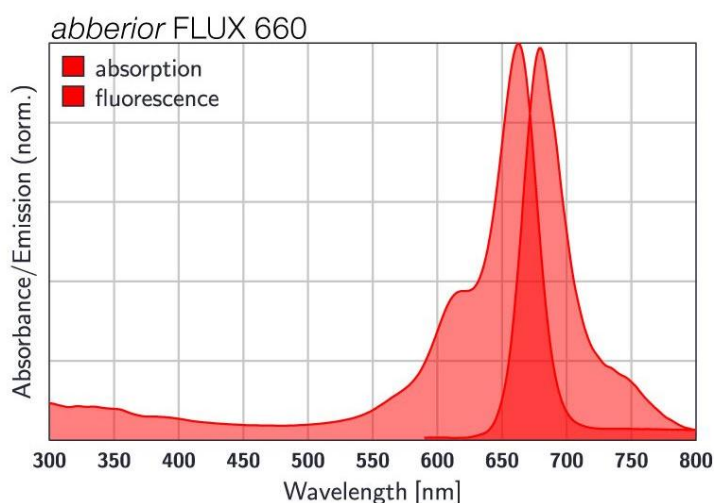

Our abberior FLUX 660 SNAP ligands can be readily used to specifically label your SNAP-tag® fusion proteins of interest within fixed cells. All our bright and photostable abberior FLUX SNAP probes are modified by a benzylguanine substrate as functional group.

This product is covered by one or more patents owned or controlled by New England Biolabs, Inc. (NEB), and is intended for Research Use Only (RUO). Commercial use of this product may require a license from New England Biolabs, Inc. For additional information or to inquire about commercial use, please contact busdev@neb.com. SNAP-tag® and CLIP-tag® are registered trademarks of New England Biolabs, Inc.

Properties

Absorption	λ_{ex} [nm]	663
Extinction Coefficient	ϵ_{max} [M ⁻¹ cm ⁻¹]	210000
Emission	λ_{em} [nm]	680
Quantum Efficiency	η_{fl} [%]	19
Fluorescence Lifetime	τ_{fl} [ns]	0.9
Correction Factor 260	CF ₂₆₀	0.13
Correction Factor 280	CF ₂₈₀	0.12
Charge	Δq	-4
Molecular Weight	MW [g/mol]	1213.3
Derivative/Conjugate	SNAP (Benzylguanine)	

Photophysical properties were measured for carboxylic acid in PBS pH 7.4.

Storage

abberior SNAP labels are shipped freeze-dried in amounts of 30 nmol per vial. Upon arrival, the substances can be stored for up to one year at -20 °C. Shortly before the staining procedure dissolve the probe in DMF or DMSO. Once dissolved the stock solutions should be kept at -20 °C, protected from light and moisture.

- Note: Depending on solvent quality the shelf-life of the stock solutions might be significantly reduced compared to the substances in their solid form – even if stored at -20 °C.