


abberior FLUX 680, goat anti-rabbit IgG, 500 µg

Item number

FX680-1002-500UG

Description

abberior FLUX 680 is a dye developed for MINFLUX imaging in the red spectral region. It combines exceptional photostability, high brightness and low background labeling with photoswitching kinetics being optimized for MINFLUX imaging. abberior FLUX 660 can be applied with excitation laser wavelengths between 630 - 690 nm together with a detection window between 650 - 750 nm. For MINFLUX imaging abberior FLUX 680 is used with MINFLUX imaging buffers such as GLOX buffer. Our abberior FLUX 680 can substitute dyes like Alexa Fluor® 680, or CF 680®. The combination of abberior FLUX 680 and abberior FLUX 640 is recommended for 2 color MINFLUX imaging.

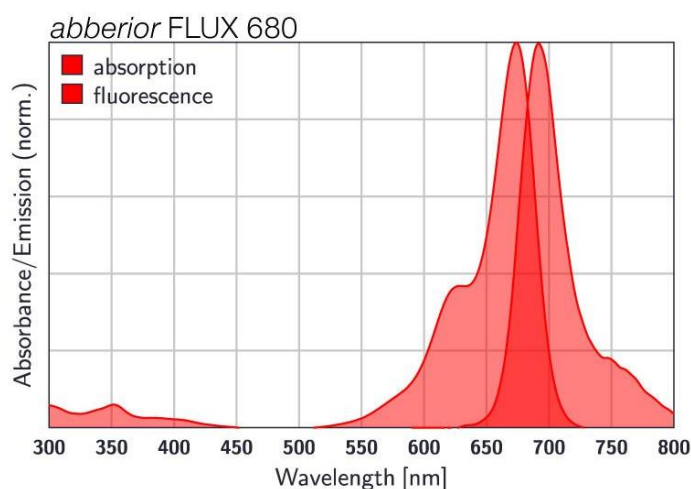

Organic fluorescent dye conjugated with polyclonal secondary anti-rabbit IgG antibody, host: goat, concentration 1 mg/ml.

Properties

Absorption	λ_{ex} [nm]	674
Extinction Coefficient	ϵ_{max} [M ⁻¹ cm ⁻¹]	250000
Emission	λ_{em} [nm]	695
Quantum Efficiency	η_q [%]	18
Fluorescence Lifetime	τ_f [ns]	0.9
Correction Factor 260	CF ₂₆₀	0.09
Correction Factor 280	CF ₂₈₀	0.1
Charge	Δq	-5
Derivative/Conjugate	anti-Rabbit IgG	

Photophysical properties were measured for carboxylic acid in PBS pH 7.4.



Hans-Adolf-Krebs-Weg 6
37077 Göttingen, Germany
E-Mail: support@abberior.com
Website: abberior.rocks

Storage

The product contains a vial of buffered solution and is shipped at room temperature. Upon arrival, the product can be stored at 4 °C for up to one month. Please split the product into smaller aliquots and store these at -20 °C to -80 °C for long-term storage of up to one year. Protect the conjugate from direct light exposure and avoid repeated freeze-thaw cycles.