


abberior STAR 488, goat anti-rabbit IgG, 500 µl (1 mg/ml)

Item number

ST488-1002-500UG

Description

abberior STAR 488 was developed for STED and confocal microscopy in the green spectral region. It can be very effectively excited with the prominent 488 nm laser line. For STED microscopy, this fluorescent dye can be efficiently used with STED laser wavelengths between 590 - 610 nm. Introducing negatively charged groups into its molecular structure making the dye excellent water soluble. Our abberior STAR 488 can substitute dyes like Oregon Green® 488, ATTO® 488 or Alexa Fluor® 488. Together with our abberior STAR 440SXP you can obtain optimal 2 color STED results with a single depletion laser. Best results are obtained with freshly prepared samples.

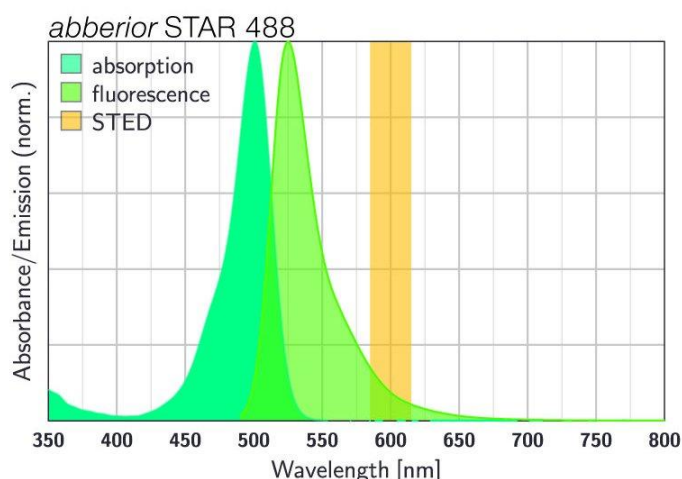

Organic fluorescent dye conjugated with polyclonal secondary anti-rabbit IgG antibody, host: goat, concentration 1 mg/ml.

Properties

Absorption	λ_{ex} [nm]	503
Extinction Coefficient	ϵ_{max} [M ⁻¹ cm ⁻¹]	65000
Emission	λ_{em} [nm]	524
Quantum Efficiency	η_f [%]	89
STED min.	$\lambda_{STED min}$ [nm]	590
STED max.	$\lambda_{STED max}$ [nm]	610
Fluorescence Lifetime	τ_f [ns]	3.9
Correction Factor 260	CF ₂₆₀	0.28
Correction Factor 280	CF ₂₈₀	0.14
Charge	Δq	-2
Derivative/Conjugate	anti-Rabbit IgG	

Photophysical properties were measured for carboxylic acid in PBS pH 7.4. Fluorescence Lifetime was measured in MeOH + 0.05% TFA.



Hans-Adolf-Krebs-Weg 6
37077 Göttingen, Germany
E-Mail: support@abberior.com
Website: abberior.rocks

Storage

The product contains a vial of buffered solution and is shipped at room temperature. Upon arrival, the product can be stored at 4 °C for up to one month. Please split the product into smaller aliquots and store these at -20 °C to -80 °C for long-term storage of up to one year. Protect the conjugate from direct light exposure and avoid repeated freeze-thaw cycles.