

abberior STAR 580, goat anti-guinea pig IgG, 500 µl (1 mg/ml)


Item number

ST580-1006-500UG

Description

abberior STAR 580 was developed for STED and confocal microscopy in the orange spectral region. Characteristic features of this dye are exceptional photostability and high fluorescent quantum yield of 90%. abberior STAR 580 can be very effectively excited with excitation wavelength between 550 - 600 nm. For STED microscopy, abberior STAR 580 can be used most efficiently with a STED laser wavelength between 750 - 800 nm. Our abberior STAR 580 can substitute dyes like Cy3, ATTO® 590 or Alexa Fluor® 568. abberior STAR 580 is the ideal partner for abberior STAR RED or abberior STAR 635P to obtain optimal 2 color STED results. Together with our abberior STAR ORANGE you get the ideal pair for FLIM experiments. Best results are obtained with freshly prepared samples.

Organic fluorescent dye conjugated with polyclonal secondary anti-guinea pig IgG antibody, host: goat, concentration 1 mg/ml.

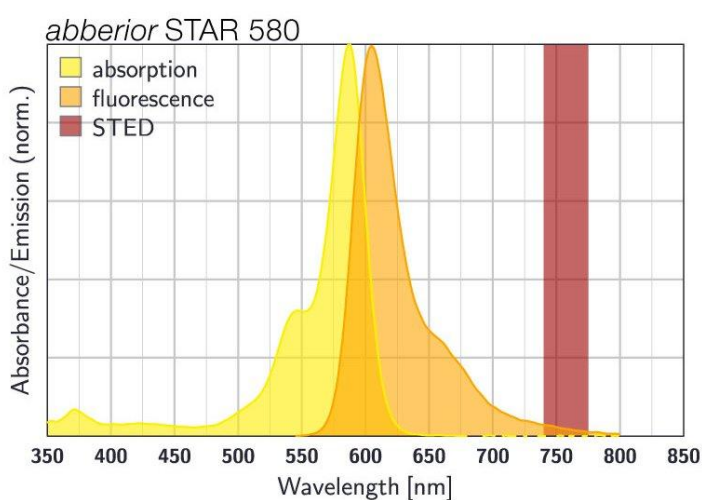

Properties

Absorption	λ_{ex} [nm]	587
Extinction Coefficient	ϵ_{max} [M ⁻¹ cm ⁻¹]	85000
Emission	λ_{em} [nm]	607
Quantum Efficiency	η_{fl} [%]	90
STED min.	$\lambda_{\text{STED min}}$ [nm]	750
STED max.	$\lambda_{\text{STED max}}$ [nm]	800
Fluorescence Lifetime	τ_{fl} [ns]	3.5
Correction Factor 260	CF ₂₆₀	0.17
Correction Factor 280	CF ₂₈₀	0.17
Charge	Δq	1
Derivative/Conjugate	anti-Guinea pig IgG	

Photophysical properties were measured for carboxylic acid in PBS pH 7.4.

Storage

The product contains a vial of buffered solution and is shipped at room temperature. Upon arrival, the product can be stored at 4 °C for up to one month. Please split the product into smaller aliquots and store these at -20 °C to -80 °C for long-term storage of up to one year. Protect the conjugate from direct light exposure and avoid repeated freeze-thaw cycles.