


abberior STAR GREEN, goat anti-guinea pig IgG, 500 µl (1 mg/ml)

Item number

STGREEN-1006-500UG

Description

abberior STAR GREEN is a novel fluorescent dye belonging to the group of fluorescein dyes. Characteristic features of this dyes are bright green fluorescent with a fluorescent quantum yield of over 80% and exceptional photostability under STED and confocal imaging conditions. abberior STAR GREEN can be very effectively excited with the 488 nm laser line. For STED microscopy, this dye can be efficiently used with a STED laser wavelength between 590 – 610 nm. Our abberior STAR GREEN is the ideal partner for abberior STAR 440L to obtain optimal 2 color STED results with a single depletion laser. abberior STAR GREEN can substitute dyes like Oregon Green® 488, ATTO® 488 or Alexa Fluor® 488. Best results are obtained with freshly prepared samples.

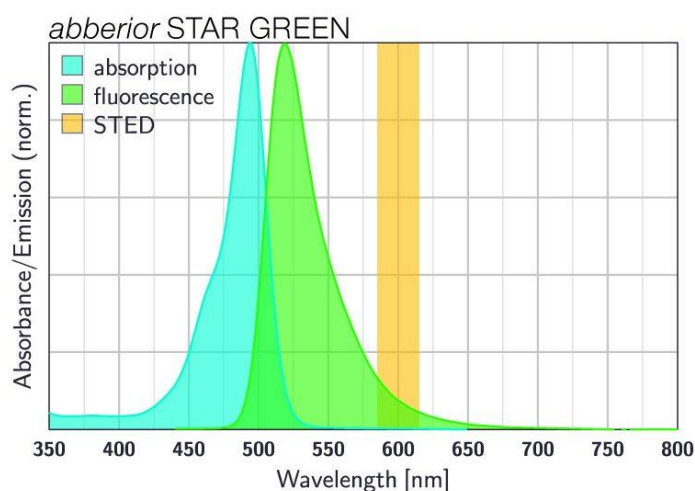

Organic fluorescent dye conjugated with polyclonal secondary anti-guinea pig IgG antibody, host: goat, concentration 1 mg/ml.

Properties

Absorption	λ_{ex} [nm]	493
Extinction Coefficient	ϵ_{max} [M ⁻¹ cm ⁻¹]	75000
Emission	λ_{em} [nm]	519
Quantum Efficiency	η_f [%]	84
STED min.	$\lambda_{STED min}$ [nm]	590
STED max.	$\lambda_{STED max}$ [nm]	610
Fluorescence Lifetime	τ_f [ns]	4.1
Correction Factor 260	CF ₂₆₀	0.48
Correction Factor 280	CF ₂₈₀	0.34
Charge	Δq	0
Derivative/Conjugate	anti-Guinea pig IgG	

Photophysical properties were measured for carboxylic acid in PBS pH 7.4.



Hans-Adolf-Krebs-Weg 6
37077 Göttingen, Germany
E-Mail: support@abberior.com
Website: abberior.rocks

Storage

The product contains a vial of buffered solution and is shipped at room temperature. Upon arrival, the product can be stored at 4 °C for up to one month. Please split the product into smaller aliquots and store these at -20 °C to -80 °C for long-term storage of up to one year. Protect the conjugate from direct light exposure and avoid repeated freeze-thaw cycles.